



**Meltis Cabinet BL261**

## At a glance :

- Furniture outline with an excellent definition
- Built-in table top
- Large storage capacity
- Optimised accessibility
- Personalization of laminated colours
- Easy maintenance, cleaning and disinfection.

## The cabinet that brings contrast to your room

The bedside cabinet BL261 has been designed to increase safety and autonomy of the residents. By using different materials and colours the furniture outline can be easily identified even for residents with impaired vision.

Both caregivers and residents can take advantage of its variety of features in terms of storage capability and options of configuration. It has a built-in table top that can be configured at right or left side of the cabinet to match with any room layout.

**Hill-Rom**

Enhancing Outcomes for Patients and Their Caregivers...

# Meltis Cabinet BL261

BL261 has a large functional storage capacity with several accessories to fulfil the needs of care institutions. The use of selected materials and removable parts are intended for an easy maintenance, clean and disinfection of the cabinet.



### Options

- Table at right side
- Table at left side
- 2 accessory bars.
- D75-mm casters, including 2 with brakes.

### Accessories

- Small basket, large basket,
- magazine holder,
- telephone holder.



### Technical specifications :


Bedside locker made of coloured melamine panels with anodised aluminium corners.

The top of the locker is covered with a protective layer made of light-beige, thermoformed ABS-PC.

The top section of the locker contains a drawer fitted with a removable plastic tray. The bottom section of the cabinet contains a large drawer fitted with two removable plastic trays.

The cast moulded polyurethane foam base, offering optimal wall protection, is fitted with 4 pivoting casters (D50 mm).

Width .....	432 mm
(with optional accessory bar) .....	446 mm
Depth .....	447 mm
Height	
(with caster - dia. Ø 50 mm) .....	907 mm
(with caster - dia. Ø 75 mm) .....	929 mm

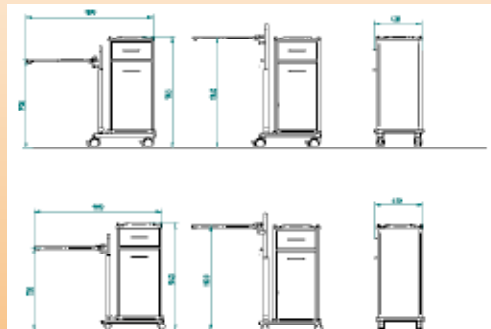

### Certification Qualité et Environnementale du site de développement et de production :

ISO 9001 (2000) LNE n° 01391 - ISO 13485 (2003) G.MED n°0851 - ISO 14001 (1996) ECOPASS n° F-563601



Enhancing Outcomes for Patients and Their Caregivers...

France.....	+33 (0)2 97 50 92 12	Suisse/Schweiz.....	+41 (0)21/706 21 30	Nordic Countries.....	+46 (0)8 564 353 60
United Kingdom.....	+44 (0)1530 41 1000	(deutschsprachig).....	.....021/706 21 38	Export.....	Contactez votre distributeur Hill-Rom local
Deutschland.....	+49 (0)211/16450-0	Österreich.....	+43 (0)2243 / 28550	.....	.....ou votre manager régional via le site web
Nederland.....	+31 (0) 347 32 35 32	Ireland.....	+353 (0)1 413 6005	.....	.....ou composez le +1 812 934 8173
Italia.....	+39 02-950541	Iberia.....	+34 (0)93 685 60 00		

Manufactured by Provita Medical GmbH. Distributed by Hill-Rom Inc.  
 Hill-Rom se réserve le droit de changer tout élément, spécification ou modèle sans notification préalable. La seule garantie dont répond Hill-Rom est la garantie expresse écrite accompagnant le contrat de vente ou de location de ses produits. © 2008 Hill-Rom Services, Inc. TOUS DROITS RÉSERVÉS.

