



Meltis Cabinet CH284

At a glance :

- Completely double-sided furniture
- A built-in tilted table top
- Perfect adaptability to any resident room
- Large storage capacity
- Ergonomic design and optimised accessibility
- Personalization of laminated colours and handles

Double-side for, double comfort

The bedside cabinet CH284 combines the comfort of an over the bed table with the large storage capacity of the Meltis cabinets.

It simplifies the work of the caregiver and it enhances the comfort and autonomy of the resident.

Meltis Cabinet CH284



- Four wheels D50 mm
- Rounded corners and edges. No sharp angles to avoid cuts and injuries.
- 0.8 mm laminated surfaces to assure perfect clean and disinfection.
- Personalization of laminated colours (24) and handles (20).
- Soft handles for visual deficient or dementia residents.
- Steel sliders in all drawers facilitate effortless and silent opening and closing.

Technical specifications :

Bedside cabinet CH284 has an integrated table top. Box and front are made of laminated particle board with ABS edging. Top of the cabinet is made of Medium-density fiberboard (MDF).

It is a double-sided bedside cabinet. It can be placed at both sides of the bed thanks to its slide-trough drawers. The base is fitted with four pivoting casters (D 50mm).

Weight 45 Kg
 Maximum load for the table top 7 Kg
 Length 586 mm
 Depth 504 mm
 Height (with D50mm casters) 870 mm
 Max. Height of table top 1120 mm
 Min. Height of table top 750 mm

Top drawer :


Length x depth x height 365x350x85 mm

Bottom drawer :

Length x depth x height 420x380x384 mm

Options :

- Bed pan holder
- Bottle holder

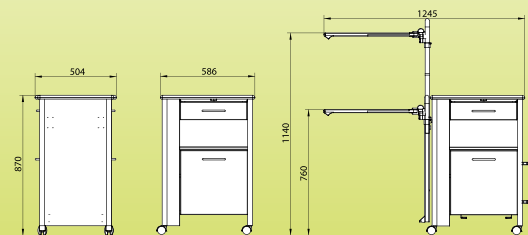

Certification Qualité et Environnementale du site de développement et de production :

ISO 9001 (2000) LNE n° 01391 - ISO 13485 (2003) G.MED n°0851 - ISO 14001 (1996) ECOPASS n° F-563601



Enhancing Outcomes for Patients and Their Caregivers...

France.....+33 (0)2 97 50 92 12
 United Kingdom.....+44 (0)1530 41 1000
 Deutschland.....+49 (0)211/16450-0
 Nederland.....+31 (0) 347 32 35 32
 Italia.....+39 02-950541

Suisse/Schweiz.....+41 (0)21/706 21 30
 (deutschsprachig).....021/706 21 38
 Österreich.....+43 (0)2243 / 28550
 Ireland.....+353 (0)1 413 6005
 Iberia.....+34 (0)93 685 60 00

Nordic Countries.....+46 (0)8 564 353 60
 Export.....Contactez votre distributeur Hill-Rom local
ou votre manager régional via le site web
ou composez le +1 812 934 8173



Manufactured by Provita Medical GmbH. Distributed by Hill-Rom Inc.
 Hill-Rom se réserve le droit de changer tout élément, spécification ou modèle sans notification préalable. La seule garantie dont répond Hill-Rom est la garantie expresse écrite accompagnant le contrat de vente ou de location de ses produits. © 2008 Hill-Rom Services, Inc. TOUS DROITS RÉSERVÉS.